**NORMANN**

User Manual / Руководство По Эксплуатации / Manual De Instrucciones / Bedienungsanleitung / Návod K Použití / Használati Utasítás / Қолдану Нұсқаулығы / Rokasgrāmata / Instrukcijų Knygėlė / Manual De Instrucțiuni

**Model: ASB-463**

GB Bathroom Scale  
RU Весы Напольные  
ES Báscula de baño  
DE Personenwaage  
CZ Osobní Váhy

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
**IMPORTANT SAFEGUARDS**

- Please read all instruction manuals before use and save it for future references.
- For home use only. Do not use for industrial or commercial purposes.
- Do not shock the scale with sudden loading or striking.
- Do not overload the scale.
- Slippery when wet.
- Do not attempt to repair, adjust or replace parts in the appliance. Repair the malfunctioning appliance in the nearest service center.
- If the product has been exposed to temperatures below 0°C for some time it should be kept at room temperature for at least 2 hours before turning it on.
  - The manufacturer reserves the right to introduce minor changes into the product design without prior notice, unless such changes influence significantly the product safety, performance, and functions.

**INSERTING THE BATTERIES**

Open the battery cover on the back of the measuring platform. Insert the Batteries into the battery case.

**Note:** If you do not intend to use this unit for a long period of time, it is advisable to remove the battery before storage.

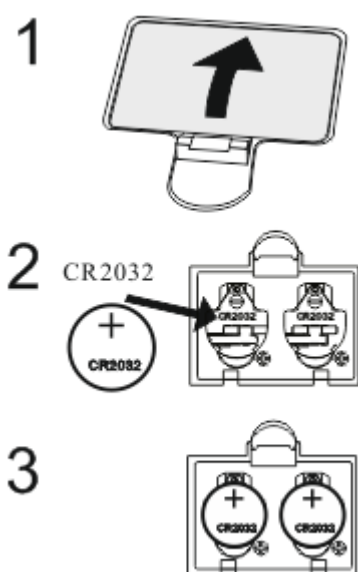**A POSITIONING THE MONITOR**

Place the measuring platform on a hard, flat surface where there is minimal vibration to ensure safe and accurate measurement.

**Note:** To avoid possible injury, do not step on the edge of the platform.

**HANDLING TIPS**

- This monitor is a precision instrument utilizing state-of-the-art technology. To keep the unit in the best condition, follow these instructions carefully.
- Do not attempt to disassemble the measuring platform.
  - Place the unit in an area free from direct sunlight, heating equipment, high humidity, or extreme temperature change.
  - Never submerge in water. Use alcohol to clean the electrodes and glass cleaner (applied to cloth first) to keep them shiny, soaps is not allowed.
  - Do not step on the platform when wet.

**USEFUL FEATURES**

- This scale allow you to measure weight, body fat percentage, water percentage, muscle percentage, bone mass and calories analysis simultaneously and easily, simply by stepping on a scale.
- Personal data (Female/Male, Age, Height) can be pre-set and stored in personal data memories.
- Automatic Shut-down function, during programming, if you do not touch any of the keys or buttons within 30 seconds. After you have completed the measuring process, the power shuts off automatically within 15 seconds.
- With Automatic Zero function.
- With Low battery and overload indication function.

**HOW TO DETERMINE BODY FAT PRECENTAGE**

**Setting and Storing Data in Memory**

The unit can be operated if data has been programmed into one of the personal data memories.

**1. Turn on the Power**

Step on the platform/Touch the SET button to turn on the unit, the unit will confirm activation; the personal data number last used will be appeared alternately.

**2. Change of unit**

Long touch the SET button, the Kg will be flash, touch the UP/DOWN button change into Lb or St

**3. Select a Personal data Code**

Firstly Press the SET button, the Personal data Code will be flashed. Press the UP/DOWN buttons to select a Personal data code. Once you reach the Personal data code you wish to use, press the SET button to confirm.

**4. Select Female or Male**

Use the UP/DOWN buttons to scroll through female, male settings, and then press the SET button to confirm

**5. Specify the Height**

Use the UP/DOWN buttons to specify Height and then press the SET button to confirm (range of user height is between 100cm-220cm).

**6. Set Age**

Use the UP/DOWN buttons to specify age. When you reach your age, press the SET button to confirm (range of user age is between 10-80 year).

**Note:** If you make a mistake before you have finished programming it, repeat the steps from Select a Personal data Code

**TAKING WEIGHT AND BODY FAT READINGS**

After programming your personal data, you are ready to take a reading.

**Select Personal Data Code**

Press the SET button turn on. Use the Up/Down buttons to select your Personal data code. The display will show the programmed data. Until the display will be show "0.0" , now step onto the platform.

**Note:** To ensure accuracy reading should be with bare-foot to step on the platform and be sure that your heels are correctly aligned with the electrodes on the measuring platform. If you step onto the platform before "0.0" appears, you will not obtain a reading. Furthermore, if you do not step onto the measuring platform within about 10 seconds after "0.0" appears, the power is shut-off automatically.

### Get Your Reading

Your weight will be shown first. Continue to stand on the platform. "0000" will appear on the display and disappear one by one from left to right. The display will then flash your body fat percentage, water percentage, bone mass, muscle percentage, calories analysis, weight alternately twice.

**Note:** Do not step off until body fat% is shown, even if you step off the platform after all the results are displayed. The alternating display of body weight and body fat percentage etc will be repeated twice, the unit will then shut-off.

### REPLACING THE BATTERIES

When the batteries run low, the Lo message will appear on the display window. When the message appears, be sure to replace the batteries immediately. Since weak batteries will affect the accuracy of your measurements. Change all the batteries at the same time with new batteries.

### ERROR MESSAGE

The following error messages may appear:

- Lo: The battery in the scales must be replaced.
- Err : The maximum load capacity is exceeded, (range of weight with 3-180kg)
- Err 2 : The body-fat and body-water percentages cannot be calculated.(If body fat percentage is less than 5% or more than 60%, the reading cannot be accurate.)

### MAINTENANCE AND CARE

- Wipe the scale with a soft cloth and cleaner, and then dry it.
- Do not use any organic solvents, aggressive chemicals and abrasives.

### STORAGE

- Complete all requirements of chapter MAINTENANCE AND CARE.
- Do not put anything on bathroom scale when not in use.
- Keep the scale in a dry cool place.

### SPECIFICATIONS

Capacity 180kg (396lb/28.8st)  
 Division 100g/0.2lb  
 LCD Display with backlight  
 Body fat percentage  
 Water percentage  
 Muscle percentage  
 Bone mass, calories  
 Automatic Zero Resetting  
 Automatic Switch off  
 Overload indication  
 Low Battery indication

### BODY FAT RATIO CHART

AGE	FEMALE(%)				MALE(%)			
	Low	Healthy	High	Obese	Low	Healthy	High	Obese
10-39	14-20	21-25	26-31	32-38	11-17	18-23	24-29	30-36
40-55	15-21	22-26	27-32	33-38	12-18	19-24	25-30	31-37
56-80	16-22	23-27	28-33	34-38	13-19	20-25	26-31	32-38

**Note:** Due to the natural differences between people this chart should be used for reference purpose only.

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